




# Isolation Gowns

Comfort drives compliance

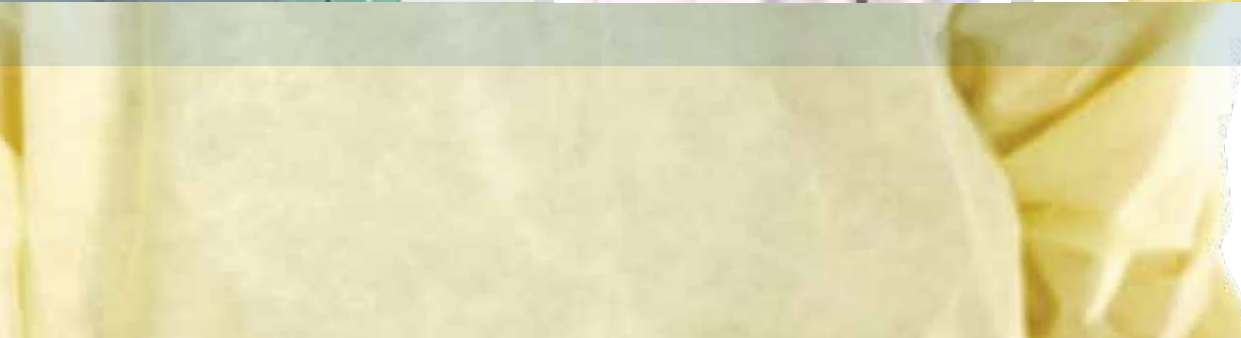

# Medline Isolation Gowns

## Comfort drives compliance

### Find your comfort zone

What makes you comfortable? Reliable protection, breathable fabric, smart features or a wide selection? When you choose Medline isolation gowns, you get all this and more.

### Protection you can trust ... in any situation

Different tasks call for different levels of fluid protection. Medline offers gowns made from several different fabrics to cover every need. Our multi-ply gowns resist fluids under low to moderate pressure, while our AAMI PB70-compliant gowns provide 360 degrees of total protection. Administering chemotherapy drugs? Choose our breathable chemo-tested gown on page 11.

### Comfort drives compliance

Who wants to wear a gown that's hot and uncomfortable? Simple logic says healthcare workers are more likely to wear appropriate protective apparel if it feels good.

### Convenient features, plenty of choices

- Tie, tape tab, over-the-head, and hook-and-loop neck closures
- Elastic wrists, thumb loops and knit cuffs
- Generous sizing, long ties
- Multiple colors: yellow, blue, green and white




### Medline's isolation gowns offer three wrist options and four neck closure options



Elastic cuffs



Thumbs Up<sup>®</sup> thumb loop

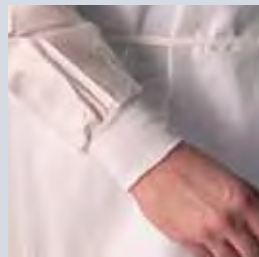

Knit cuffs



Hook-and-loop neck closure for easy adjustment



Tape tab neck closure for easy dressing and removal

We created easy-to-use icons to help identify location of use based on our recommendations.


|   |   |   |  |
|---|---|---|--|
| Decontamination      | Emergency  | Food Services  | Airborne Precautions  |
| Intensive Care Unit  | Nursing    | Visitors       | Contact Precautions   |
| Maternity            | Lab        | Radiology      | Fluid Resistance      |

# AAMI Gowns

## AAMI Level 1

### Lightweight Polypropylene Isolation Gowns

8

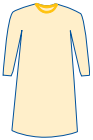


#### Yellow, Elastic Wrists

|            |               |       |
|------------|---------------|-------|
| NONLV100   | Regular/Large | 50/cs |
| NONLV100XL | X-Large       | 50/cs |

## AAMI Level 2

### Medium weight Polypropylene Isolation Gowns

8

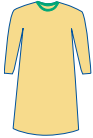

#### Yellow, Elastic Wrists

|             |               |        |
|-------------|---------------|--------|
| NONLV200M   | Medium        | 100/cs |
| NONLV200    | Regular/Large | 100/cs |
| NONLV200XL  | X-Large       | 100/cs |
| NONLV200XXL | XX-Large      | 100/cs |

#### Yellow, Thumb Loop Wrists, Over-the-Head Neck

|            |               |        |
|------------|---------------|--------|
| NONLV240   | Regular/Large | 100/cs |
| NONLV240XL | X-Large       | 100/cs |

#### Yellow, Elastic Wrists, Individually Folded

|             |               |       |
|-------------|---------------|-------|
| NONLV200F   | Regular/Large | 40/cs |
| NONLV200FXL | X-Large       | 40/cs |

#### Yellow, Thumb Loop/Individually Folded

|              |               |       |
|--------------|---------------|-------|
| NONLV200TL   | Regular/Large | 40/cs |
| NONLV200TLXL | X-Large       | 40/cs |

## AAMI Level 3

### Kazzmere Premium Gowns

9




#### Blue, Knit Cuffs

|            |               |       |
|------------|---------------|-------|
| NONLV325   | Regular/Large | 50/cs |
| NONLV325XL | X-Large       | 50/cs |

### Heavyweight Multi-Ply Premium Gowns

9

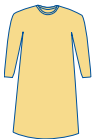

#### Yellow, Elastic Wrists

|            |               |        |
|------------|---------------|--------|
| NONLV315   | Regular/Large | 100/cs |
| NONLV315XL | X-Large       | 100/cs |

### Microporous Breathable Laminate Gowns

10



#### White, Knit Cuffs

|            |               |       |
|------------|---------------|-------|
| NONLV350   | Regular/Large | 50/cs |
| NONLV350XL | X-Large       | 50/cs |

# Non-AAMI Gowns

## Multi-Ply Standard Isolation Gowns

4



#### Yellow, Medium weight, Elastic Wrists

|             |               |        |
|-------------|---------------|--------|
| NON27SMS2   | Regular/Large | 100/cs |
| NON27SMS2XL | X-Large       | 100/cs |

#### Blue, Medium weight, Elastic Wrists

|             |               |        |
|-------------|---------------|--------|
| NON27SMS3   | Regular/Large | 100/cs |
| NON27SMS3XL | X-Large       | 100/cs |

#### Yellow, Lightweight, Elastic Wrists

|            |               |       |
|------------|---------------|-------|
| NON27236   | Regular/Large | 50/cs |
| NON27236XL | X-Large       | 50/cs |


#### Yellow, Medium weight Level 2 Material, Thumb Loop Wrists\*

|             |               |        |
|-------------|---------------|--------|
| NON27550    | Regular/Large | 100/cs |
| NON27550XL  | X-Large       | 100/cs |
| NON27SMS5   | Regular/Large | 100/cs |
| NON27SMS5XL | X-Large       | 100/cs |

\*Made with same material as used in level 2 gowns. Not tested against AAMI standards.

## Spunbond Polypropylene Cover Gowns

5

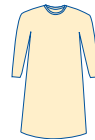

#### Yellow, Elastic Wrists

|         |               |       |
|---------|---------------|-------|
| CRI4000 | Regular/Large | 50/cs |
| CRI4001 | X-Large       | 50/cs |
| CRI4003 | XXX-Large     | 50/cs |

#### Yellow, Over-the-Head, Thumb Loop Wrists

|         |               |       |
|---------|---------------|-------|
| CRI4200 | Regular/Large | 50/cs |
| CRI4201 | X-Large       | 50/cs |

## Coated Polypropylene Gowns

6



#### Blue, Knit Cuffs

|             |               |       |
|-------------|---------------|-------|
| NON27116    | Regular/Large | 50/cs |
| NON27116XL  | X-Large       | 50/cs |
| NON27116XXL | XX-Large      | 50/cs |

#### Blue, Elastic Wrists

|            |               |       |
|------------|---------------|-------|
| NON27114   | Regular/Large | 50/cs |
| NON27114XL | X-Large       | 50/cs |

#### Blue, Over-the-Head/Knit Cuffs

|             |               |       |
|-------------|---------------|-------|
| NON27279OH  | Regular/Large | 50/cs |
| NON27279OHL | X-Large       | 50/cs |

#### Blue, Over-the-Head/Thumb Loop

|               |               |       |
|---------------|---------------|-------|
| NON27279OHTL  | Regular/Large | 50/cs |
| NON27279OHTLX | X-Large       | 50/cs |

#### Yellow, Knit Cuffs

|           |               |       |
|-----------|---------------|-------|
| CRI4010CF | Regular/Large | 50/cs |
| CRI4011CF | X-Large       | 50/cs |

#### Yellow, Elastic Wrists

|         |               |       |
|---------|---------------|-------|
| CRI4010 | Regular/Large | 50/cs |
| CRI4011 | X-Large       | 50/cs |

## Polyethylene (Thumb Loop Style) Gowns

7



|          |               |       |
|----------|---------------|-------|
| CRI5000  | Regular/Large | 75/cs |
| CRI5010  | Universal     | 75/cs |
| CRI5001  | X-Large       | 75/cs |
| CRI5003  | XXX-Large     | 75/cs |
| NONTH150 | Regular/Large | 75/cs |
| NONTH180 | Universal     | 75/cs |
| NONTH200 | X-Large       | 75/cs |

\*Above gowns are individually folded

## Prevention Plus Gowns

11



#### Blue, Knit Cuffs

|             |               |       |
|-------------|---------------|-------|
| NON27457    | Regular/Large | 30/cs |
| NON27457XL  | X-Large       | 30/cs |
| NON27457XXL | XX-Large      | 30/cs |

KEY: Unless noted individually folded, all gowns are flat packed in bags of 10.

# Multi-Ply Isolation Gowns

More protection than standard cover gowns

## Medium weight Multi-Ply Fabric



Multi-Ply  
Fabric  
Swatch  
Above and  
Below

## Lightweight Multi-Ply Fabric

### Multi-ply material

The material is made from densely packed meltblown inner layers sandwiched between strong, spunbond outer layers. These gowns are available in two weights: medium weight and lightweight.

### Gown features

Gowns feature neck and waist ties with elastic wrists.

### Fluid-resistant protection

Both lightweight and medium weight provide fluid-resistant protection.

Medium Weight Gown

Lightweight Gown

### Low-linting

Because these gowns are made from polypropylene, they are naturally low in lint.



| Medium Weight Multi-Ply Polypropylene Protection and Comfort Levels |      |        |      |
|---|------|--------|------|
|   | Good | Better | Best |
| Protection  |      |        |      |
| Comfort   |      |        |      |



| Lightweight Multi-Ply Polypropylene Protection and Comfort Levels |      |        |      |
|---|------|--------|------|
|   | Good | Better | Best |
| Protection  |      |        |      |
| Comfort   |      |        |      |



# Polypropylene Cover Gowns

For basic protection



## Polypropylene Fabric Swatch

### Polypropylene cover gowns

Gowns made from polypropylene offer basic cover protection.

### Gown features

Gowns are available with either neck ties and elastic wrists or apron-style, over-the-head neck and thumb loop wrists.

### Fluid-repelling protection



Gowns are able to repel light amounts of fluids.

### Low-linting

Because these gowns are made from polypropylene, they are naturally low in lint.

**Note:** These gowns are constructed of a fabric that offers basic cover protection and does not meet the AAMI PB70 Standard. For gowns that meet AAMI standards, please refer to pages 8-10.



#### Polypropylene Protection and Comfort Levels

|            | Good  | Better | Best |
|------------|---|--------|------|
| Protection | [Progress bar showing approximately 25% completion] |        |      |
| Comfort    | [Progress bar showing approximately 75% completion] |        |      |



non-AAMI

# Coated Isolation Gowns

High fluid protection and comfort

*Medium weight  
Coated  
Fabric*



Coated  
Fabric  
Swatch  
Above and  
Below

*Heavyweight  
Coated  
Fabric*

## Coated material

Made from polypropylene and polyethylene coated together. Choose from two weights: medium weight and heavyweight.

## Three options for wrists

These gowns feature comfortable latex-free knit or elastic wrists with neck and waist ties. Also available with thumb loop wrists.

## Protection from fluids

These gowns provide great fluid protection.

## Low-linting

Polypropylene fabric makes these gowns naturally low in lint.

**Note:** *These gowns are constructed of a fabric that offers high fluid protection, yet due to seam construction, the gowns do not meet the AAMI PB70 Standard. For gowns that meet AAMI standards, please refer to pages 8-10.*



Coated Protection and Comfort Levels

|            | Good   | Better | Best |
|------------|--|--------|------|
| Protection |  |        |      |
| Comfort    |   |        |      |



# Multi-Ply AAMI PB70 Isolation Gowns

Multiple layers offer higher levels of protection

**Medium weight  
Multi-Ply  
Fabric**



**Multi-Ply  
Fabric  
Swatch  
Above and  
Below**

**Lightweight  
Multi-Ply  
Fabric**

## Multiple layers of protection

Our technology utilizes multiple layers of fabric to provide a high level of fluid repellency. The fabric features densely packed meltblown layers sandwiched between strong, spunbond outer layers. Gowns are available in two weights for two levels of protection.

## Easy to wear

The tape tab neck closure and over-the-head neck style makes these gowns easy to don and remove.

## Meets AAMI Levels 1 & 2

The gowns meet the Association for the Advancement of Medical Instrumentation (AAMI) PB70 standard, "Liquid barrier performance and classification of protective apparel and drapes intended for use in health care facilities."

*Medium Weight Gown – Level 2*

*Lightweight Gown – Level 1*

## Color-coded for AAMI protection level

The neckline on these gowns are color-coded, to correspond to the AAMI barrier performance level.

*Green Neck Binding – Level 2 Gown*

*Yellow Neck Binding – Level 1 Gown*

## Low-linting

Because these gowns are made from polypropylene, they are naturally low in lint.



|            | Good | Better | Best |
|------------|------|--------|------|
| Protection |      |        |      |
| Comfort    |      |        |      |



|            | Good | Better | Best |
|------------|------|--------|------|
| Protection |      |        |      |
| Comfort    |      |        |      |

LEVEL 1

LEVEL 2



# Premium AAMI PB70 Isolation Gowns

Strong and comfortable

**Kazzmere  
Fabric**



**Heavyweight  
Multi-Ply**

## Level 3 Isolation Gowns

Breathable, comfortable gowns that offer a high level of fluid protection.

### Convenient Features

Different options available for donning and style:  
*Kazzmere* – latex-free knit cuffs, hook-and-loop neck closure

*Heavyweight Multi-Ply* – elastic wrists, tape tab neck closure

### Protection - AAMI Approved Level 3



The gowns meet the Association for the Advancement of Medical Instrumentation (AAMI) PB70 standard, “Liquid barrier performance and classification of protective apparel and drapes intended for use in health care facilities” - **Level 3**

### Color coded for AAMI protective level

The neckline on these gown is color coded purple to correspond to the AAMI barrier performance level – **Level 3**

### Low-Linting

These isolation gowns are naturally low in lint.



Kazzmere Protection and Comfort levels

|            | Good                               | Better | Best |
|------------|------------------------------------|--------|------|
| Protection | [Progress bar from Good to Better] |        |      |
| Comfort    | [Progress bar from Good to Better] |        |      |



Heavyweight Multi-Ply Protection and Comfort levels

|            | Good                               | Better | Best |
|------------|------------------------------------|--------|------|
| Protection | [Progress bar from Good to Better] |        |      |
| Comfort    | [Progress bar from Good to Better] |        |      |



LEVEL 3



# Prevention Plus Gowns

Highly impervious protection; breathable microporous technology



## Prevention Plus Fabric Swatch

### Advanced technology delivers protection, breathability and comfort

Prevention Plus gowns offer a high level of protection while providing increased breathability and comfort. The gown is composed of three layers: a patent-pending hybrid film laminated between two layers of bi-component spunbond material. Because moisture vapor is thousands of times smaller than the smallest known virus, the microporous film allows moisture vapor to escape while preventing transmission of fluids. This effect increases comfort without sacrificing protection.



Prevention gown's breathable, microporous film technology protects against the visually undetectable while still allowing moisture vapor to escape

### Sonically-sealed seams prevent strike-through

Prevention gowns feature sonically-sealed seams to prevent strike-through where the sleeves are joined to the body of the gown.



Sonically-sealed seams provide a higher level of protection

### Material meets standard test methods

Impervious breathable material passed ASTM 1670 and ASTM 1671 testing methods.

### Chemotherapy drug tested

Prevention Plus gowns are tested against the top chemotherapy drugs. See chart below.

The Prevention Plus Material has been tested against the following chemotherapy drugs per ASTM F739-99a:

|                            |                     |
|----------------------------|---------------------|
| Dacarbazine                | Etoposide (Toposar) |
| Carmustine (BiCNU)         | Paclitaxel (Taxol)  |
| Cyclophosphamide (Cytoxan) | Thio-Tepa           |
| Doxorubicin Hydrochloride  | Mitomycin C         |
| 5-Fluorouracil             | Vincristine Sulfate |
| Cisplatin                  |                     |

Test data available upon request



#### Prevention Plus Protection and Comfort Levels

|            | Good                                  | Better | Best |
|------------|---------------------------------------|--------|------|
| Protection | [Progress bar showing 'Better' level] |        |      |
| Comfort    | [Progress bar showing 'Good' level]   |        |      |



LIT171



LIT821

Ask your Medline representative or call 1-800-MEDLINE for help selecting the right protective apparel for your facility.

Request our *Protective Apparel brochure* to view our complete protective apparel line, including generously sized footwear, protective headwear, gowns that meet AAMI PB70 requirements, protective multi-purpose gowns, protective coveralls, lab coats and lab jackets, *LIT 171*.

To find out more about the wide variety of Over-the-Head isolation gowns that we offer, ask for a copy of our Over-the-Head Personal Protection Gowns brochure, *LIT 821*.



**Medline Industries, Inc.**

One Medline Place  
Mundelein, IL 60060

**Medline United States**  
**1-800-MEDLINE (633-5463)**

www.medline.com | info@medline.com

**Medline Canada**  
**1-800-396-6996**

www.medline.ca | canada@medline.com

**Medline México**  
**01-800-831-0898**

www.medlinemexico.com | mexico@medline.com