



## AREACARE RANGE


A comprehensive range of static foam pressure relieving surfaces designed to maintain skin integrity and help prevent pressure injury, while delivering maximum comfort.

# AreaCare Range

As specialists in the supply of pressure care and manual transfer equipment solutions designed to help prevent pressure wound and skin integrity issues, Novis support a range of clinical healthcare products designed for hospitals, aged care facilities, rehabilitation providers and the home care community.

With a long history delivering effective pressure care management solutions, we are committed to delivering successful patient outcomes through product innovation, professional collaboration and clinical excellence.


Our product portfolio has been engineered with the goal of enhancing the comfort and safety of patients; reducing the time and concerns of care staff; and improving clinical outcomes for care facilities.

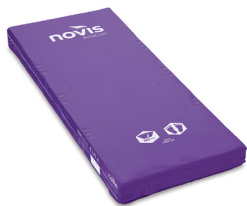

## AreaCare Pro

- Multi-density, three layer, zoned design for maximum comfort
- Hinged side walls which aid bed profiling and help prevent falls
- Castellated memory foam top to aid pressure redistribution and articulates with movement

<b>SIZE</b>	Single
<b>CAPACITY</b>	220 kg

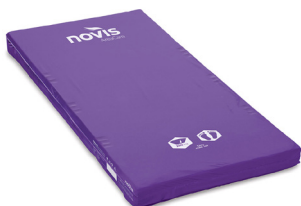

## AreaCare Plus

- Multi-density, three layer, zoned design for maximum comfort
- Hinged side walls which aid bed profiling and help prevent falls
- Extra height (18cm) King Single version for greater weight capacity and support

<b>SIZE</b>	Single	King Single
<b>CAPACITY</b>	220 kg	250 kg



## AreaCare

- Multi-density, three layer design
- Memory foam top layer for optimal pressure redistribution
- Foot end zipper with waterfall skirt

<b>SIZE</b>	Single	King Single
<b>CAPACITY</b>	200 kg	220 kg



## Underlay

- Provides “bottoming out” protection for alternating pressure overlays
- Minimises height of overall support system

<b>SIZE</b>	Single
<b>CAPACITY</b>	NA*



## Extension

- 200mm extension for beds with sleeping platform extension
- Memory foam top layer provides additional pressure relief under ‘at risk’ heel zone

<b>SIZE</b>	Single or King Single (18 cm)
<b>CAPACITY</b>	NA*

# An extensive range of features to offer every patient optimum care and comfort

Multi-density foam layers for maximum strength and comfort – made from premium material, each foam layer is made of varying densities to maximise the redistribution of pressure.

## Top Cover

Latex free, multi-stretch polyurethane coated top cover for reduced shear and friction. Vapour permeable, bacteriostatically treated, fire retardant and water resistant. Double-layered waterfall skirt provides protection from moisture ingress.

Cover includes an asset management panel and turning schedule (where applicable).

## Infection Control

Top cover seams are welded and/or tape sealed for maximum infection control.

AreaCare Pro and AreaCare Plus feature a 360° double-layered waterfall skirt to direct fluids away from mattress base.

## Hinged Base

Six hinge system to conform to profiling beds.

## Side Bolster

Safety side bolster provide a firm surface for reduced risk of falls when sitting on edge of mattress. Includes a six hinge system to conform to profiling beds.

## Top Layer

High density, temperature sensitive, memory foam conforms to the body's shape; working in combination with the higher resistance support layer to disperse pressure evenly. Special immersion, slow recovery foam foot zone for the vulnerable heel area.

## Inner Cotton Cover

Foam covered by a removable inner cotton cover for added protection.

## Support Layer

High resilience foam for maximum pressure redistribution.

## Core Layer

Firm core layer and edges provide form, strength and stability for maximum safety and reduce risk of falls.

## Base Cover

Hard-wearing, welded, and water resistant base cover material. AreaCare Pro and AreaCare Plus feature a 360° YKK zipper for easy removal and two carrying handles per side for easy handling.



Find out more about the AreaCare range at <https://novis.com.au/category/brands/areacare/>



# NOVIS®

enabling independence



**AreaCare Pro**

**AreaCare Plus**

**AreaCare Plus King**

**AreaCare**

**AreaCare King**

**Extension**

**Extension King**

<b>Length</b>	1960 mm	1960 mm	1960 mm	1960 mm	1960 mm	200 mm	200 mm
<b>Width</b>	865 mm	865 mm	1050 mm	865 mm	1050 mm	865 mm	1050 mm
<b>Height</b>	150 mm	150 mm	180 mm	150 mm	150 mm	150 mm	180 mm
<b>Layers</b>	3	3	3	3	3	2	2
<b>Side walls</b>	✓	✓	✓	✗	✓	✗	✗
<b>360° Zipper</b>	✓	✓	✓	✗	✗	✗	✗
<b>SWL</b>	220 kg	220 kg	250 kg	200 kg	220 kg	NA	NA
<b>Infection Control Seams</b>	✓	✓	✓	✓	✓	✓	✓
<b>Hinged</b>	✓	✓	✓	✗	✓	✗	✗
<b>Castellation/Convolution</b>	✓	✗	✗	✗	✗	✗	✗
<b>Waterfall Skirt</b>	✓	✓	✓	✓ (foot end only)	✓ (foot end only)	✗	✗
<b>Warranty</b>	6 year body; 2 year cover						
<b>Code</b>	BSPS3300	BSPS3000	BSPS3001	BSPS3003	BSPS3002	BSPS5300	BSPS5301