

STERIGLIDE™ CANNULA

50%
Better Gliding

STERIGLIDE™ The cannula that glides.

STERIJECT™ NEEDLE

67%
Better Flow Rate

ULTRA Thin Wall needles for dermal filling.

14%
Thinner than 33G

THE INVISIBLE NEEDLE™
The thinnest needle ever made.


0.08^{ml}
BoNT Savings

33G low dead space needle.

THE INVISIBLE NEEDLE™

LOW DEAD SPACE

Sharp Needles

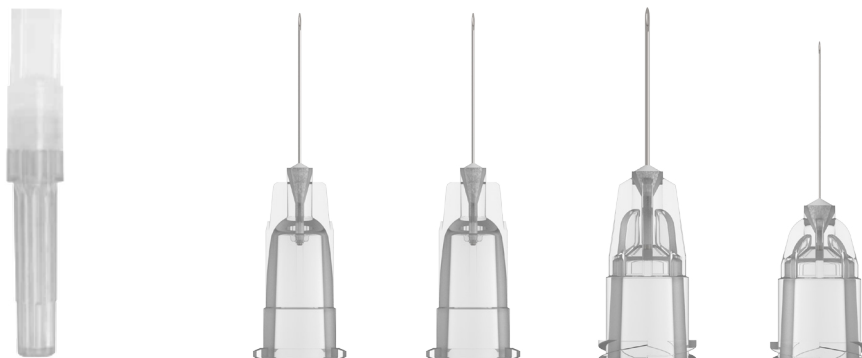

Hub Features	PRE	PRC	HPC	LDS	Hub Benefits
Hard Pack Capsule	X	X	X	X	Preventing Needle Injury
Hard Polymer		X	X	X	No Flexing
Threading			X	X	No Pop-Off Risk
Low Dead Space				X	No Product Loss
Needle Up to 14% thinner than 33G				X	Minimal Patient Trauma
Recommended Use	BoNT + Meso	Filler	Filler	BoNT	

THE INViSIBLE NEEDLE™.


THE INViSIBLE NEEDLE™ is the first needle of its kind. Technology advancements have made it possible to make this needle 33% thinner than a 30G needle (used for insulin syringes) and 14% thinner than the TSK 33G needle. This extremely small diameter of THE INViSIBLE NEEDLE™ results in a mere painless experience for patients and is ideal for BoNT injections.



NEW

Low dead space hub.

The **Low Dead Space needle hub** is designed to reduce the effective space left behind in the needle hub as close to zero. This new patented LDS hub design results in savings of up to 0.08ml which leads to significant cost savings when injecting high cost drugs, such as botulinum toxin.

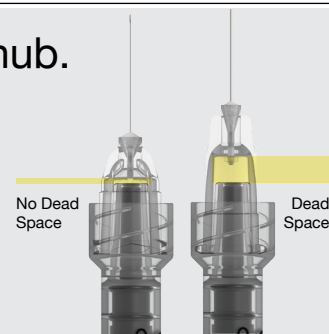

BoNT Needles / LDS Hub

Code	Description	Inch	UTW*	Box
LDS-02009	THE INViSIBLE NEEDLE x 9mm**	3/8"		100
LDS-33013	Low Dead Space Hub 33G x 13mm**	1/2"		100
LDS-33009	Low Dead Space Hub 33G x 9mm**	3/8"		100

BoNT & Mesotherapy Needles / PRE Hub

Code	Description	Inch	UTW*	Box
PRE-33013	Regular Hub 33G x 13mm	1/2"		100
PRE-33004	Regular Hub 33G x 4mm	3/16"		100
PRE-32013	Regular Hub 32G x 13mm	1/2"		100
PRE-32009	Regular Hub 32G x 9mm	3/8"		100
PRE-32006	Regular Hub 32G x 6mm	1/4"		100
PRE-32004	Regular Hub 32G x 4mm	3/16"		100
PRE-31013	Regular Hub 31G x 13mm	1/2"		100
PRE-30006I	Regular Hub 30G x 6mm	1/4"	✓	100
PRE-30004I	Regular Hub 30G x 4mm	3/16"	✓	100

Filler Needles / PRC Hub

Code	Description	Inch	UTW*	Box
PRC-30038I	Control Hub 30G x 38mm	1 1/2"	✓	100
PRC-30025I	Control Hub 30G x 25mm	1"	✓	100
PRC-30013I	Control Hub 30G x 13mm	1/2"	✓	100
PRC-27038I	Control Hub 27G x 38mm	1 1/2"	✓	100
PRC-27019I	Control Hub 27G x 19mm	3/4"	✓	100
PRC-27013I	Control Hub 27G x 13mm	1/2"	✓	100
PRC-25013I	Control Hub 25G x 13mm	1/2"	✓	100
PRC-22050I	Control Hub 22G x 50mm	2"	✓	100


Filler Needles / HPC Hub

Code	Description	Inch	UTW*	Box
HPC-30013I	Advanced Hub 30G x 13mm	1/2"	✓	320
HPC-27013I	Advanced Hub 27G x 13mm	1/2"	✓	320
HPC-23025I	Advanced Hub 23G x 25mm	1"	✓	320

* UTW = ULTRA Thin Wall ** For exclusive use with Soft-Ject Luer Lock syringe

All other sizes on request

Cannulas

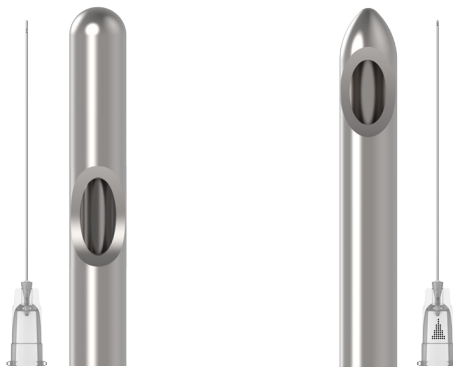

Features	CSH	STERiGLIDE™	Benefits
Safety PRC Hub	X	X	No Pop-Off Risk
ULTRA Thin Wall	X	X	Best Product Flow Rate
Dome Shaped Tip		X	Easy Cannula Introduction
Nearest To Tip Filler Delivery		X	More Accurate Placement
Proprietary Surface Treatment		X	Better Gliding
Clear Side-port Marking		X	Better Handling


STERiGLIDE™ cannula.

The biggest advantage of cannulas over needles is the reduced risk of bruising during filler injections. This advantage has resulted in a rapid adoption of the cannula technique worldwide, making cannulas the product of choice for selective procedures, e.g. the tear through region.



Other advantages of cannulas include the ability to treat a larger area with only one injection point, better comfort, less trauma, and more acceptability by the patient.

After TSK's first cannula introduction with the CSH Cannula, TSK introduced their award winning STERiGLIDE cannula in 2012. The STERiGLIDE outperforms any other cannula available and remains to lead the market as the golden standard.

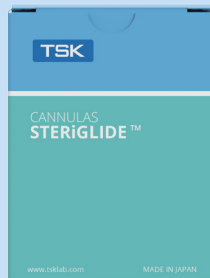


Proven 50% better gliding.

When cannulas were first introduced the biggest challenges were to penetrate the cannula into the incision entry hole, to move it easily under the skin and to inject the filler accurately.

The design and function of a regular cannula lacked these benefits, thus a new cannula needed to be developed.

TSK successfully introduced the STERiGLIDE™ to the market and overcame these challenges. The STERiGLIDE™ features a newly designed tip, a proprietary coating, and a clear side hole marking. These new improvements allowed for up to 50% better gliding, the ability to control the direction of the cannula as well as the ability to place the filler more precisely.

■ 22G, STERiGLIDE Cannula
■ 22G, Regular CSH Cannula
■ 22G, Sharp needle

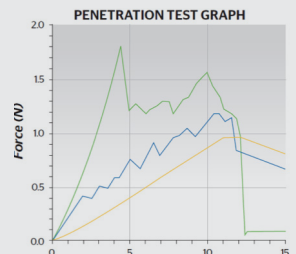

Cannulas / STERiGLIDE™


Code	Description	Inch	UTW*	Box
PRC-30025ISG	STERiGLIDE Cannula 30G x 25mm	1"	✓	20
PRC-27050ISG	STERiGLIDE Cannula 27G x 50mm	2"	✓	20
PRC-27038ISG	STERiGLIDE Cannula 27G x 38mm	1½"	✓	20
PRC-27025ISG	STERiGLIDE Cannula 27G x 25mm	1"	✓	20
PRC-25050ISG	STERiGLIDE Cannula 25G x 50mm	2"	✓	20
PRC-25038ISG	STERiGLIDE Cannula 25G x 38mm	1½"	✓	20
PRC-22070ISG	STERiGLIDE Cannula 22G x 70mm	2¾"	✓	20
PRC-22050ISG	STERiGLIDE Cannula 22G x 50mm	2"	✓	20

Cannulas / CSH

Code	Description	Inch	UTW*	Box
PRC-30025ICSH	CSH Cannula 30G x 25mm	1"	✓	50
PRC-27050ICSH	CSH Cannula 27G x 50mm	2"	✓	50
PRC-27038ICSH	CSH Cannula 27G x 38mm	1½"	✓	50
PRC-27025ICSH	CSH Cannula 27G x 25mm	1"	✓	50
PRC-25050ICSH	CSH Cannula 25G x 50mm	2"	✓	50
PRC-25038ICSH	CSH Cannula 25G x 38mm	1½"	✓	50
PRC-22070ICSH	CSH Cannula 22G x 70mm	2¾"	✓	50
PRC-22050ICSH	CSH Cannula 22G x 50mm	2"	✓	50

*UTW = ULTRA Thin Wall

Syringes & Accessories

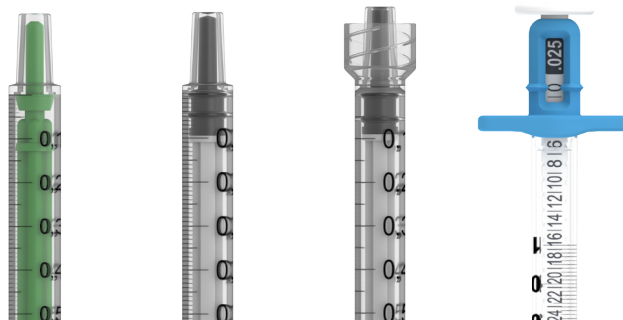

Features	A1	S1	L1	3Dose	Benefits
Low Dead Space	X	X	X	X	No Product Loss
Rubber Stopper		X	X	X	Smooth Flow
Luer Lock			X	X	No Pop-Off Risk
Dose Clicker System				X	Save up to 20% BoNT
Unit Indicator				X	Inject per unit
Tactile, Audible				X	100% Accurate Dose Injection

3Dose: 100% accurate injection.

The two most critical variables during BoNT injections are the amount and the anatomical place where the BoNT is being injected.

The new innovative BoNT injection system, the 3Dose, features an easy click system that accurately controls the injection amount preventing over/under dosage.


With a unique clicker system it provides audible and tactile feedback on the precise amount of BoNT being injected. This eliminates the need to focus on how much to inject, thereby allowing the focus to solely be on where to inject.



Prevent needle injury.

We continue to look for accessories that make our product more convenient, safer and smarter.

The safety aspect of needle puncture is always present and for this reason a needle holder is offered where our needle capsule packaging fits in tightly and securely.

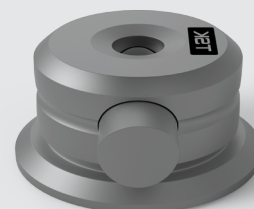

Syringes

Code	Description	ml	Box
DOSE3-33013	3Dose System + LDS 33G x 13mm	1	20
HSW-L1	3-part Soft-Ject Luer Lock 1ml	1	100
HSW-S1	3-part Soft-Ject Luer Slip 1ml	1	100
HSW-A1	2-part Norm-Ject Luer Slip 1ml	1	100

Accessories

Code	Description	Box
TNH-1	TSK Needle Holder	1